
# SIDE TO SIDE WITH MASSAGE

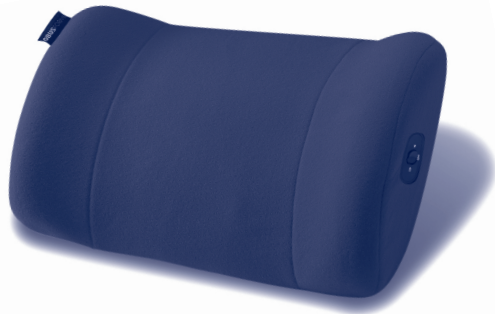English, Français

## CONGRATULATIONS

on your purchase of an ObusForme Side to Side with Massage. The ObusForme Side to Side with Massage gently supports the natural curvature of your lower back, preventing the tendency to slouch and helps to enhance overall sitting posture. By encouraging proper posture the ObusForme Side to Side with Massage helps to minimize fatigue by reducing the strain and energy demand on your body.

The ObusForme Side to Side with Massage can provide relief from pain aggravated by poor posture including back pain, neck pain, shoulder tension and headaches. By reducing spinal compression it can also ease the pain caused by spinal disc problems and sciatica.

## HOW TO USE YOUR OBUSFORME SIDE TO SIDE WITH MASSAGE

Small, lightweight and portable, your ObusForme Side to Side with Massage is ideal for use anywhere you sit – at home, in the office and on the go. Simply place the Side to Side with Massage against the back of your chair or seat in an upright position, wherever it feels the most comfortable for you. Use the attached adjustable elastic strap to secure it in place with ease. To activate the massage option, adjust the switch to (I) low intensity or (II) high intensity depending on your preferred setting

## USE & CARE TIPS FOR YOUR OBUSFORME SIDE TO SIDE WITH MASSAGE

Your ObusForme Side by Side with Massage is made with high-density polyurethane foam and other synthetic materials. If the surface of the ObusForme Side to Side with Massage becomes soiled, it can be wiped gently with a damp cloth, mild soap and a sponge, or a soft brush. Do not rub the product excessively or place it in a washing machine. Avoid soaking the product as this may damage the material.

### BATTERY REPLACEMENT

1. Locate the battery compartment on the back of the Side-To-Side with Massage.
2. Remove the exhausted batteries from the Side-To-Side with Massage.
3. Insert four, new AA (LR6) alkaline batteries as indicated inside the battery compartment.

### Helpful Hints:

- For longer battery life, we recommend use of alkaline batteries.
- Fit the tab on the battery compartment cover into the slot in the foam and close the cover.
- Do not mix old and new batteries.
- Do not mix different types of batteries: alkaline, standard (carbon-zinc) or rechargeable (nickel-cadmium).
- Remove the batteries during long periods of non-use. Always remove exhausted batteries from the Side To Side with Massage. Battery leakage and corrosion can cause damage. Dispose of batteries safely.
- Never short circuit the battery terminals.
- Non-rechargeable batteries are not to be recharged.
- Only batteries of the same or equivalent type as recommended in the “Battery Replacement” instructions are to be used.
- If removable rechargeable batteries are used, they are only to be charged under adult supervision.
- Rechargeable batteries are to be removed from the Side To Side with Massage before they are charged.
- If you use a battery charger, it should be examined regularly for damage to the cord, plug, enclosure and other parts. Do not use a damaged battery charger until it has been properly repaired.

## OBUS FORME LIMITED WARRANTY

ObusForme guarantees all items are free from defects in workmanship & materials for a period of time between the original purchase date & that stated below. This guarantee applies when items are used for the purpose intended. Items will be repaired/replaced (at our option), with new, refurbished parts/products &/or substitutes, if the ORIGINAL purchaser has completed & returned the Warranty Registration within 30 days of purchase & includes original receipt. Shipping, customs, duties & taxes must be PRE-PAID TO & FROM ObusForme by the PURCHASER. This warranty gives you rights that vary by province/state. This warranty may change.

### WHAT IS **NOT** COVERED

Wear & tear, aging (including foam discoloration, flattening, density, consistency), accidental damages, alterations, mishandling, faulty adjustment, misuse, improper care, power damage, rental use, discontinued items, service by anyone other than ObusForme, shipping damages, over inflation, neglect, items sold ‘as is’ or damage due to natural acts are NOT covered.

### WARRANTY TIME FRAME

**Backrest Frame: *Lifetime***  
(cover, foam, lumbar pad & other Backrest parts/materials are NOT covered)

**Seat Frame: *1 year***  
(cover, foam & other Seat parts/materials are NOT covered)

**Back Therapy: *1 year***  
(Back Therapy includes Backlife, Back Belts, Back Packs, Drivers Seats)

**Sleep/Foot/Muscle Therapy: *1 year***  
(pillow cases are NOT covered)

**Electrical Parts: *1 year***  
(this includes wires, adaptors, plugs and other electrical parts/components)

### HOW TO OBTAIN WARRANTY SERVICE

Please obtain a Return Authorization Number & instructions prior to sending your item or it may be denied. Please inform Customer Service by:

**Mail:**  
**HoMedics Group Canada**  
344 Consumers Road  
Toronto, Ontario, Canada M2J 1P8

Tel: (416) 785-1386 Fax: (416) 785-5862  
Toll Free: 1-888-225-7378  
8:30 a.m. to 5:00 p.m., Mon - Fri EST  
[www.obusforme.com](http://www.obusforme.com)

**OBUSFORME®**  
THE BODYCARE EXPERTS™

Corporate Member of ACE

ASSOCIATION OF CANADIAN ERGONOMISTS

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ObusForme is committed to providing you with optimal relief and comfort. To serve you better in the future, we would like to know if we have fulfilled our commitment. Please complete and return this Questionnaire to help us better meet your needs.

We aggregate this information and use it internally for research and marketing purposes only. We do not disclose personal information to any third parties.

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